



Disobliteration System

Reperfusion

Improving Outcomes

Better Economics

To become the **first line therapy** in treating peripheral vascular occlusions by developing a **versatile** solution aimed at rapid reperfusion without the need for surgery or thrombolytic drug therapies.

A percutaneous, one-size fits all thrombectomy/embolectomy system designed to replace surgical and drug therapies

Ideal for acute limb ischemia, organized thrombus, thrombosed grafts, embolic events, and post atherectomy debris



COMING SOON!
Large Bore
DVT Solutions

Features & Benefits

7F System

- Non surgical approach eliminates lengthy ICU time and surgical complications
- System compatibility promotes follow therapies to treat underlying issues

One size fits all (peripheral)

- System designed to treat above the knee and below

Control Sheath

- 7F Sheath with 10mm arrests flow which promotes a reduction in distal embolization
- Maintains access to treat potential underlying lesions

Disobliteration

(disö-blitě-rāshön)

Removal of a blood clot or atherosclerotic plaque from an obstructed blood vessel.

Atraumatic mechanism of action

- Balloon sweep is a proven means of clot removal and proven gentle on vessels
- Eliminates large volume blood loss associated with aspiration
- Reduces or eliminates drug therapy which potentiates a reduction of common bleeding complications associated with lytics
- Promotes a reduction in distal embolization common with aspiration and ultrasonic therapies
- Avoids tubular necrosis, acute renal insufficiency, and potential renal failure often associated with rapidly moving maceration or aspiration technologies

No need for capital equipment

- Eliminates added cost, maintenance, and eliminates connecting tubing and other accessories

Elegant catheter based design

- Easily adoptable by interventionalists
- Rx balloon nitinol core allows for tremendous pushability

Call today

954.483.4525

www.ichorvascular.com

Linked in



Product Details

Included with ICHOR Reperfusion System (ICH7FR)

1. 7F x45cm Control Sheath with 10mm Occlusion Balloon
2. 7Fx90cm Guide Catheter with 10mm Expandable Nitinol Funnel
3. 014 Rapid Exchange Catheter with 10mm Balloon and Nitinol Core
 - 035 Sheath Dilator
 - (2) Hemostatic Valves w/ side port
 - (2) 1cc Syringes


Funnel Capture Slide

- Easy of deployment

7Fx90cm Guide Catheter

with 10mm Expandable Nitinol Funnel / Capture

- One size fits most anatomy
- Tight weave nitinol retrieval for maximum material removal
- Minimal wall contact reduces potential for vessel wall damage
- Manual aspiration side port to further aid in vessel clearance

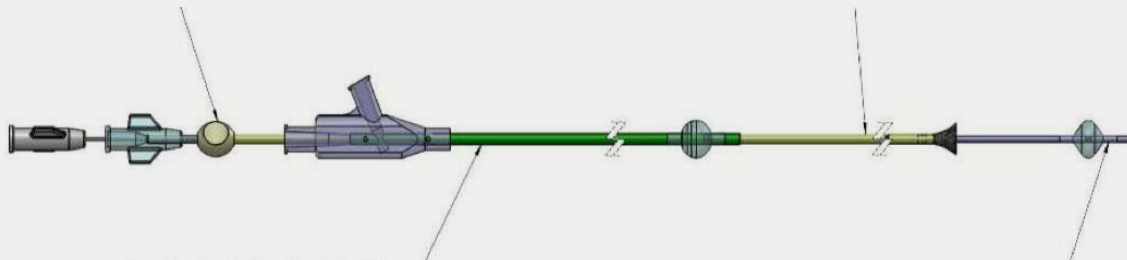

7Fx45cm Sheath with 10mm Occlusion Balloon

- Reduces potential distal embolization
- Flow reversal ideal with supplemental aspiration
- Anchors the system
- Maintains access for multiple passes
- Atraumatic

014 Rapid Exchange Catheter with 10mm Balloon and Nitinol Core

- Ideal for use with shorter wire systems
- Superior pushability
- Tapered tip for excellent trackability
- Ideal for use in wide range of vessels

Ordering Information

ICHOR Reperfusion System..... ICH7FR
 Quantity5/bx

954.483.4525

www.ichorvascular.com

Linked in